




Veraseal[®] 3

Featuring Patented AIR[°]gel[®] Technology

The Veraseal 3 Disposable mask brings AIR[°]gel[®] technology to acute care for a giant step forward in patient comfort. Three color coded elbow designs make choosing the right mask for CPAP/Bi-level, single or dual limb circuit NIV systems simple. The Veraseal 3 fits quickly and easily while enhancing patient comfort.

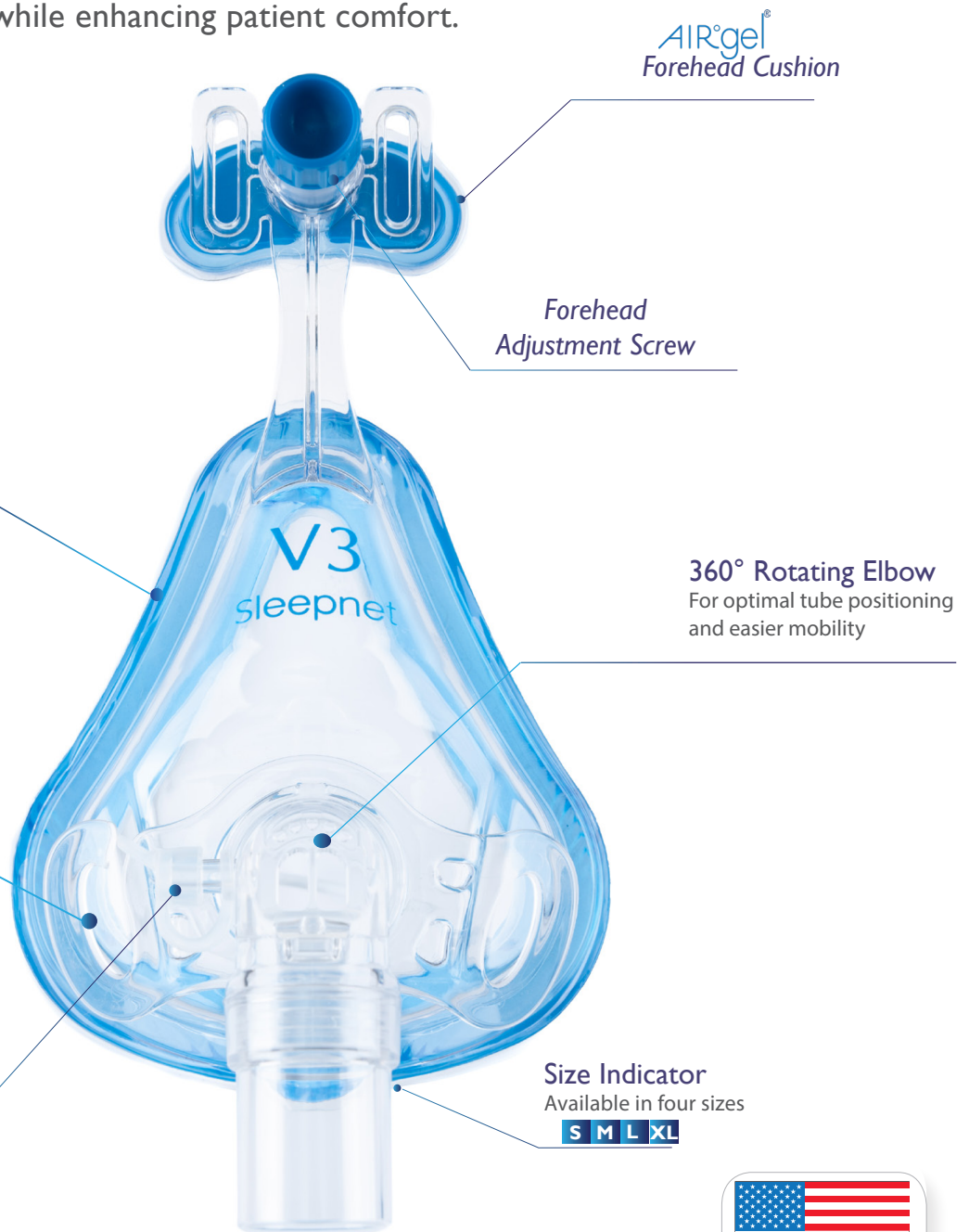

AIR[°]gel[®]
Forehead Cushion

Forehead
Adjustment Screw

AIR[°]gel[®]

Maximum comfort,
minimum weight.
Our patented
AIR[°]gel minimizes
pressure points and
skin irritation.

360° Rotating Elbow
For optimal tube positioning
and easier mobility

Quick-Release
Slot

Lightweight, breathable
headgear slides off easily
with one hand.

Oxygen Port

Allows easy delivery of supplemental
oxygen.

Size Indicator

Available in four sizes

S M L XL



What makes our Veraseal 3 Mask different?

- Designed with minimal intentional leaks and elbows that work with CPAP/Bi-level and single or dual limb circuits.
- Adjustment Screw - provides a wide, continuous range of adjustment ensuring a custom fit while removing pressure from the bridge of the nose.
- FitRight 7-Day Headgear - maintains a secure and comfortable fit during treatment.
- Available in four sizes with simple and breathable quick-release headgear.

Featuring Patented AIR^ogel[®] Technology

Our patented AIR^ogel with an ergonomic design and luxurious feel reduces pressure while maintaining an effective seal.



A full line of disposable masks to meet your NIV needs, color coded for easy identification



Vented Mask with an anti-asphyxia valve for use with CPAP or bi-level devices



Nonvented Mask with an anti-asphyxia valve for use with single limb circuits



Nonvented Mask for use with dual limb circuits

Veraseal 3 Masks with FitRight 7-Day Headgear Specifications:

Product Number	Style	Leak Port	Entrianment Valve	Size	Mask Weight	Mask Dimensions
50843	Vented	✓	✓	Small	109.2g	3.60"D x 4.10"W x 7.25"H
50844	Vented	✓	✓	Medium	114.3g	3.60"D x 4.30"W x 7.40"H
50845	Vented	✓	✓	Large	119.4g	3.60"D x 4.38"W x 7.60"H
50846	Vented	✓	✓	Extra Large	120.1g	3.60"D x 4.42"W x 7.80"H
50949	Non-Vented	✗	✗	Small	108.6g	3.70"D x 4.10"W x 7.25"H
50950	Non-Vented	✗	✗	Medium	113.7g	3.70"D x 4.30"W x 7.40"H
50951	Non-Vented	✗	✗	Large	118.8g	3.70"D x 4.38"W x 7.60"H
50952	Non-Vented	✗	✗	Extra Large	119.5g	3.70"D x 4.42"W x 7.80"H
50953	AAV Non-Vented	✗	✓	Small	109.1g	3.65"D x 4.10"W x 7.25"H
50954	AAV Non-Vented	✗	✓	Medium	114.2g	3.65"D x 4.30"W x 7.40"H
50955	AAV Non-Vented	✗	✓	Large	119.3g	3.65"D x 4.38"W x 7.60"H
50956	AAV Non-Vented	✗	✓	Extra Large	120g	3.65"D x 4.42"W x 7.80"H



Does not contain natural rubber latex

Contact us today.

General enquires email
info@sleepnetcorp.com or
 call Call 1 (603) 758 - 6600

Visit www.sleepnetmasks.com to learn more.

Sleepnet Corporation, 5 Merrill Industrial Drive,
 Hampton, New Hampshire 03842 USA